Xeva Tutorial

Arvind Mer 2017-01-17

Load Xeva

Load Xeva library and data.

```
library(Xeva)
data(lpdx)
head(modelInfo(lpdx))
```

```
##
                                    model.id
                                                                  donor
## PHLC1106 P5.501.A1.1 PHLC1106 P5.501.A1.1 11101S-213RC-312S(F)-412S-
## PHLC1106_P5.504.A4.1 PHLC1106_P5.504.A4.1 11101S-213RC-312S(F)-412S-
## PHLC1106_P5.506.B1.1 PHLC1106_P5.506.B1.1 11101S-213RC-312S(F)-412S-
## PHLC1106_P5.507.B2.1 PHLC1106_P5.507.B2.1 11101S-213RC-312S(F)-412S-
## PHLC1106_P5.508.B3.1 PHLC1106_P5.508.B3.1 11101S-213RC-312S(F)-412S-
## PHLC1106_P5.511.C1.1 PHLC1106_P5.511.C1.1 11101S-213RC-312S(F)-412S-
                             dob sex
                                         PHLC biobase.id patient.id
## PHLC1106_P5.501.A1.1 Aug31.14
                                  F PHLC1106
                                                PHLC1106
                                                           PHLC1106
## PHLC1106_P5.504.A4.1 Aug31.14
                                  F PHLC1106
                                                PHLC1106
                                                           PHLC1106
                                                           PHLC1106
## PHLC1106_P5.506.B1.1 Aug31.14
                                 F PHLC1106
                                                PHLC1106
## PHLC1106_P5.507.B2.1 Aug31.14
                                 F PHLC1106
                                                PHLC1106
                                                           PHLC1106
## PHLC1106_P5.508.B3.1 Aug31.14
                                  F PHLC1106
                                                PHLC1106
                                                           PHLC1106
## PHLC1106 P5.511.C1.1 Sep14.14
                                   F PHLC1106
                                                           PHLC1106
                                                PHLC1106
```

Models which belongs to same batch are in one list which is stored in expDesign slot. For example


```
print(batchNames(lpdx))
```

```
##
     PHLC1106 P5
                     PHLC111_P7
                                   PHLC119 P5
                                                  PHLC153_P6
                                                                 PHLC181 P7
## "PHLC1106_P5"
                   "PHLC111_P7"
                                 "PHLC119_P5"
                                                "PHLC153_P6"
                                                               "PHLC181 P7"
##
      PHLC189_P5
                     PHLC191_P5
                                   PHLC191_P7
                                                  PHLC196_P5
                                                                 PHLC215_P5
    "PHLC189 P5"
                                 "PHLC191 P7"
                                                "PHLC196 P5"
##
                   "PHLC191 P5"
                                                               "PHLC215 P5"
##
      PHLC229 P6
                    PHLC235 P4
                                   PHLC655 P7
                                                   PHLC82 P5
    "PHLC229 P6"
                   "PHLC235 P4"
                                 "PHLC655 P7"
                                                 "PHLC82 P5"
```

To calculate angle between the treatment and control samples of this batch

```
batchNames <- batchNames(lpdx)
expDesign <- expDesign(lpdx, batchNames[1])
ang <- calculateAngle(lpdx, expDesign, treatment.only = TRUE, plot=TRUE)</pre>
```



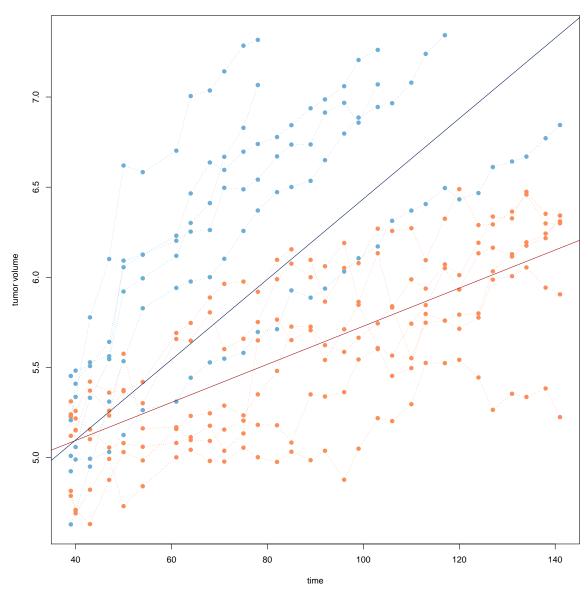


print(ang)


```
## $PHLC1106_P5
## [1] 0.671583
```


```
#par(mfrow=c(5,3))
for(I in batchNames)
{
   expDesign <- expDesign(lpdx, I)
   ang <- calculateAngle(lpdx, expDesign, treatment.only = TRUE, plot=TRUE)
# print(ang)
}</pre>
```



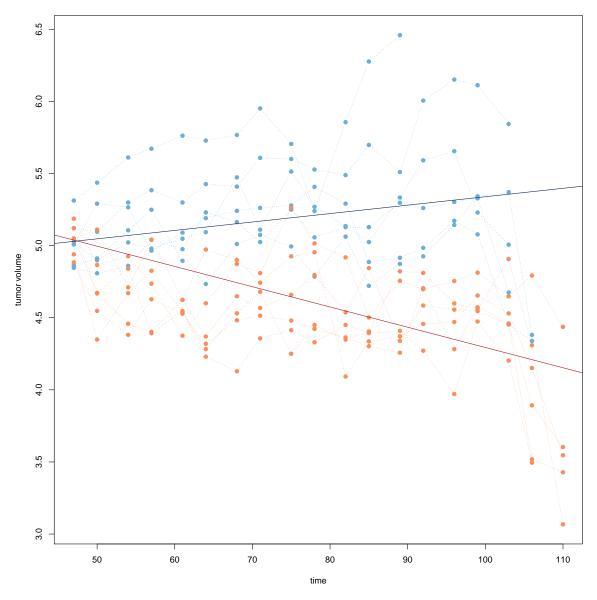




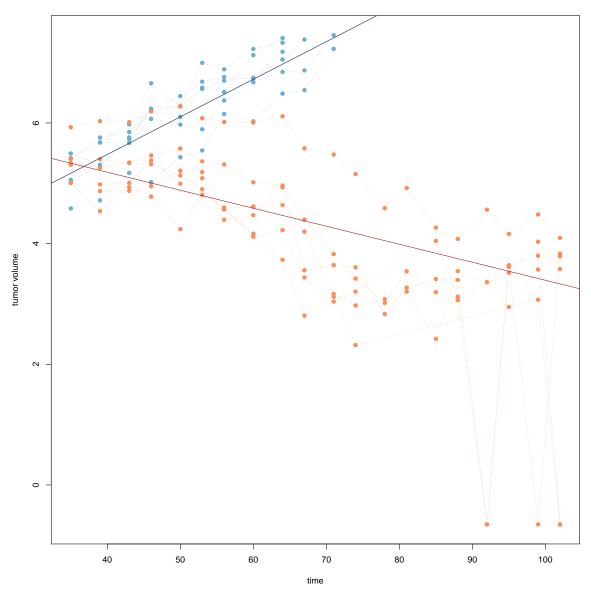


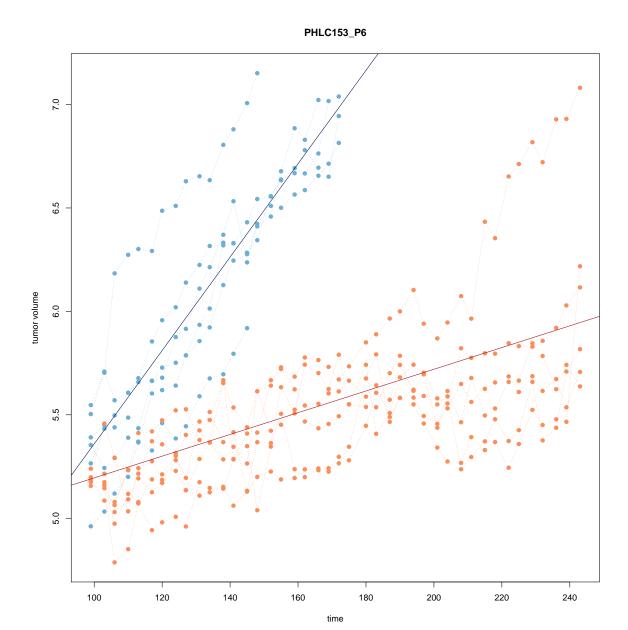




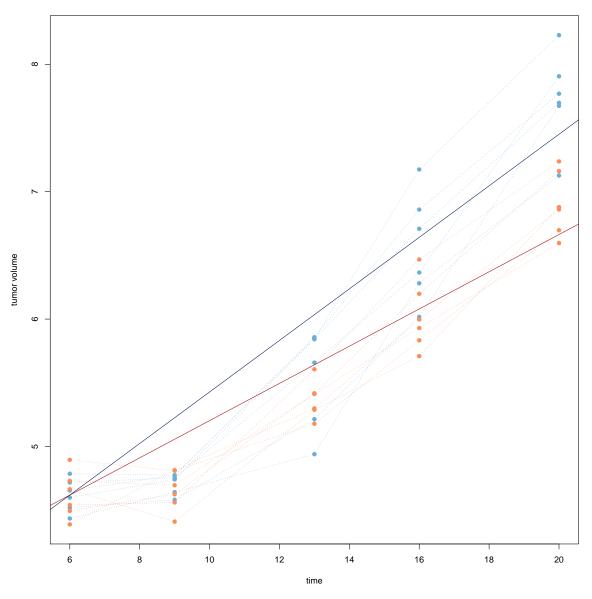


PHLC189_P5

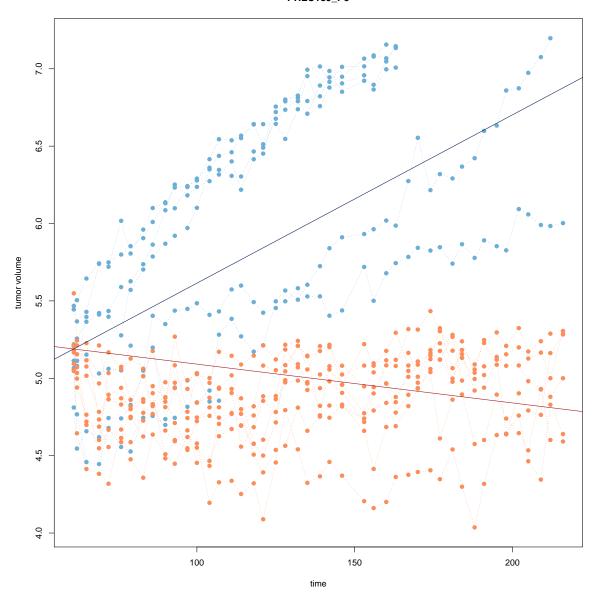


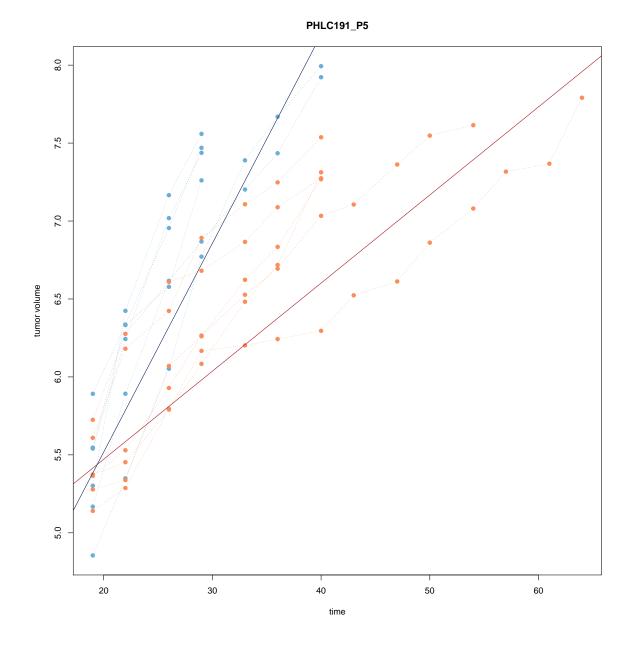




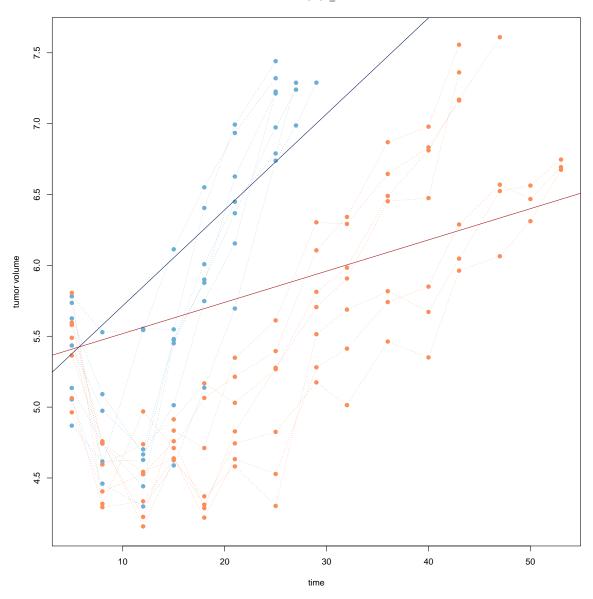




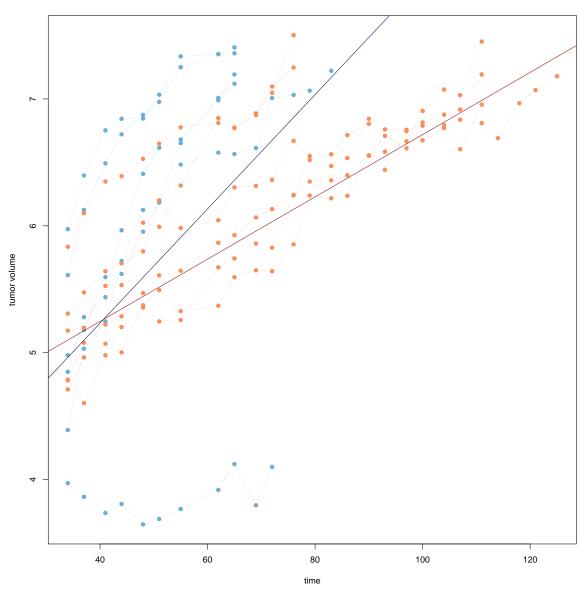





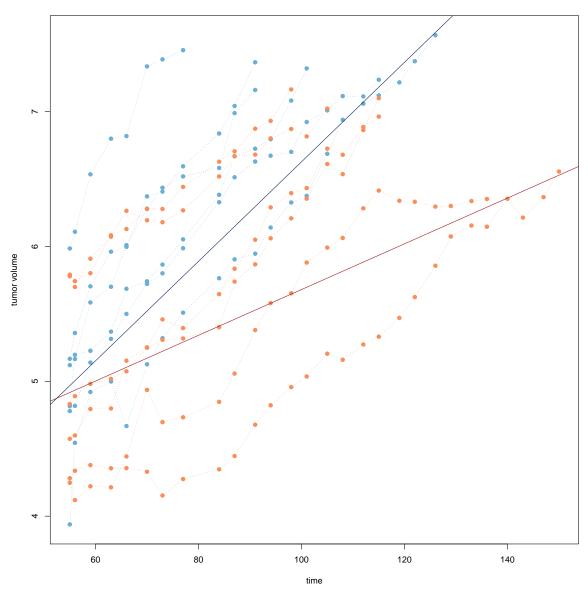

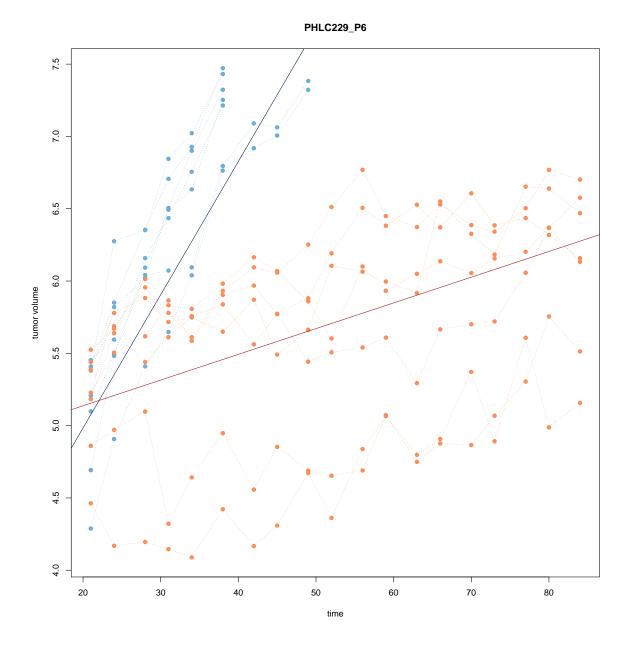



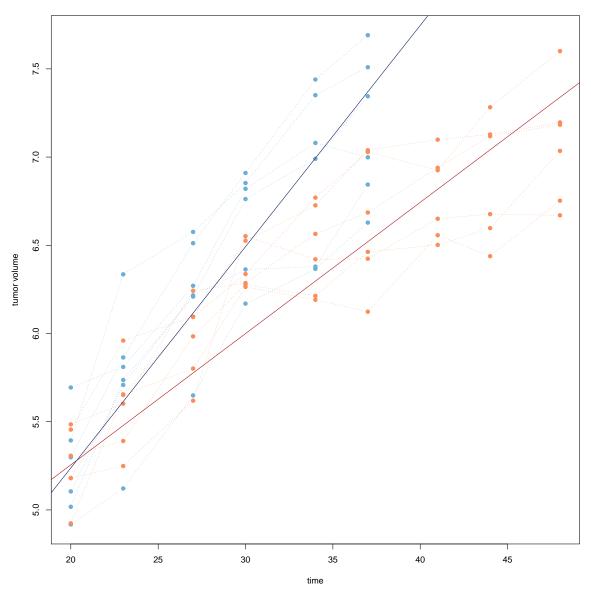

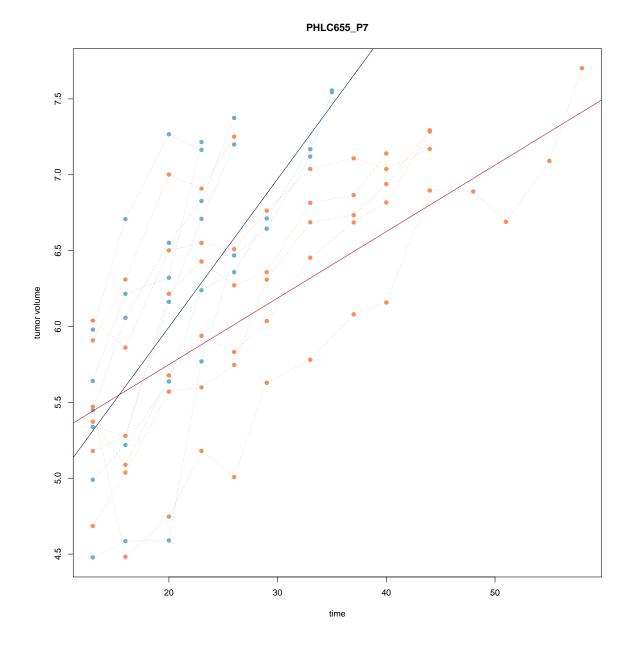

PHLC82_P5

Summarize Response of PDXs Get slop of each model and combine summarize all model slop which belongs to same patient by "mean" $^{\prime\prime}$

150 time 200

250

Get angle between treatment and control model ids. For each batch it will give one angle value

100

```
lpdx_angle <- summarizeResponse(lpdx, response.measure = "angle")</pre>
```

Get mutation expression profile

50

4.5

```
ldxe_mut <- getMolecularProfiles(lpdx, data.type="mutation")</pre>
print(ldxe_mut)
## ExpressionSet (storageMode: lockedEnvironment)
## assayData: 10483 features, 12 samples
     element names: exprs
## protocolData: none
## phenoData
     sampleNames: PHLC1106 PHLC111 ... PHLC82 (12 total)
##
     varLabels: PHLC.ID X.ID
     varMetadata: labelDescription
##
## featureData
##
     featureNames: NOC2L ISG15 ... RNF128 (10483 total)
     fvarLabels: probe.Id
##
##
    fvarMetadata: labelDescription
## experimentData: use 'experimentData(object)'
## Annotation: MUT
```

The sample names in expression set are called biobase.id in model slot. Sample names from the expression set canb be be mapped to individual PDX model.ids as

```
# get sample names
library(Biobase)
## Loading required package: BiocGenerics
## Loading required package: parallel
##
## Attaching package: 'BiocGenerics'
## The following objects are masked from 'package:parallel':
##
##
       clusterApply, clusterApplyLB, clusterCall, clusterEvalQ,
##
       clusterExport, clusterMap, parApply, parCapply, parLapply,
##
       parLapplyLB, parRapply, parSapply, parSapplyLB
## The following objects are masked from 'package:stats':
##
##
       IQR, mad, xtabs
## The following objects are masked from 'package:base':
##
##
       Filter, Find, Map, Position, Reduce, anyDuplicated, append,
```

as.data.frame, cbind, colnames, do.call, duplicated, eval, evalq, get, grep, grepl, intersect, is.unsorted, lapply,

lengths, mapply, match, mget, order, paste, pmax, pmax.int,

table, tapply, union, unique, unsplit, which, which.max,

pmin, pmin.int, rank, rbind, rownames, sapply, setdiff, sort,

##

##

##

##


##

which.min


```
## Welcome to Bioconductor
##
##
       Vignettes contain introductory material; view with
##
       'browseVignettes()'. To cite Bioconductor, see
       'citation("Biobase")', and for packages 'citation("pkgname")'.
##
sn <- Biobase::sampleNames(ldxe_mut)</pre>
smap <- mapModelSlotIds(lpdx, id=sn, id.name = "biobase.id", map.to = "model.id")</pre>
head(smap)
##
                        biobase.id
                                               model.id
## PHLC1106_P5.501.A1.1
                          PHLC1106 PHLC1106_P5.501.A1.1
## PHLC1106_P5.504.A4.1
                         PHLC1106 PHLC1106_P5.504.A4.1
## PHLC1106_P5.506.B1.1 PHLC1106 PHLC1106_P5.506.B1.1
## PHLC1106_P5.507.B2.1 PHLC1106 PHLC1106_P5.507.B2.1
## PHLC1106_P5.508.B3.1 PHLC1106 PHLC1106_P5.508.B3.1
```

PHLC1106_P5.511.C1.1 PHLC1106 PHLC1106_P5.511.C1.1

Create mRECIST plot for PDXE Lung Cancer data

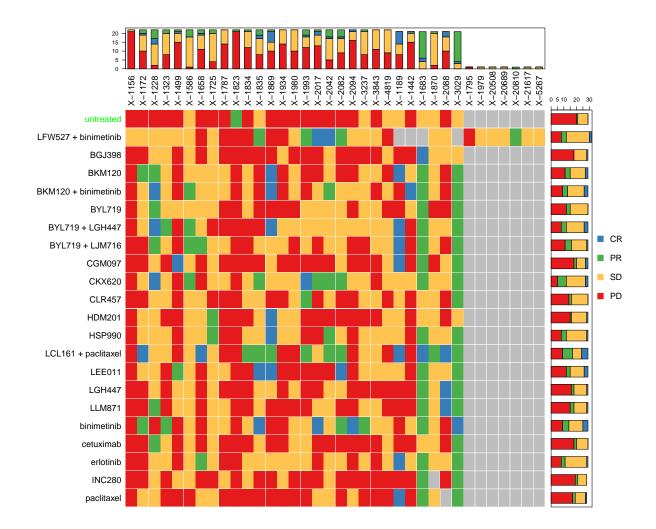

```
#pdf(file="DATA-raw/mRECIST_plot_BRCA.pdf", width=12, height=10)
brDF = df[df$tumor.type=="BRCA", ]
plotmRECIST(brDF, groupBy = "biobase.id", control.name = "untreated")
```

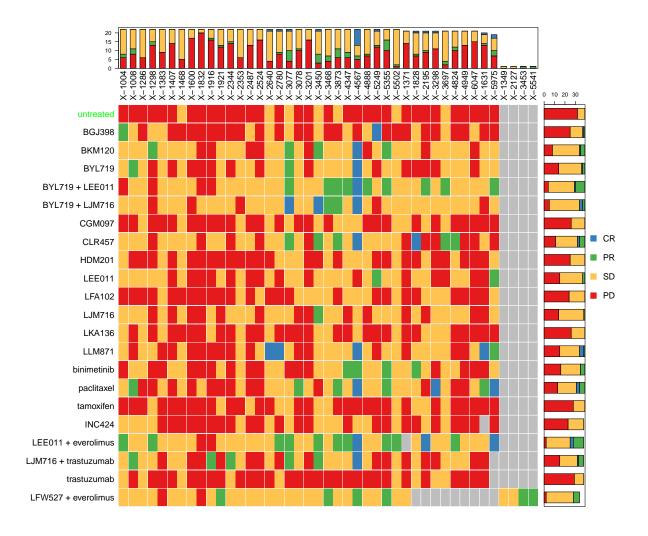

#dev.off()


Creat mR vs slop bar-plot

```
slope=c(); mR=c()
for(dn in rownames(lung_pdxe_slope))
{
    for(pi in colnames(lung_pdxe_slope))
    {
        v = c(lung_pdxe_slope[dn,pi], lung_pdxe_mR[dn,pi])
        if(!is.na(v[1]) & !is.na(v[1]))
        { slope = c(slope,v[1]); mR=c(mR,v[2]) }
    }
}

df = data.frame(mR= mR, slope= as.numeric(slope), stringsAsFactors = FALSE)
df$mR= factor(df$mR, c("CR", "PR", "SD", "PD"))

colPalette = c("#377eb8", "#4daf4a", "#fec44f", "#e41a1c")
#pdf(file="DATA-raw/boxplot_lungCancer.pdf", width=12, height=10)
boxplot(slope~mR, data=df, col=colPalette,
    main="mRECIST vs slope", xlab="mRECIST", ylab="slope")
```

mRECIST vs slope

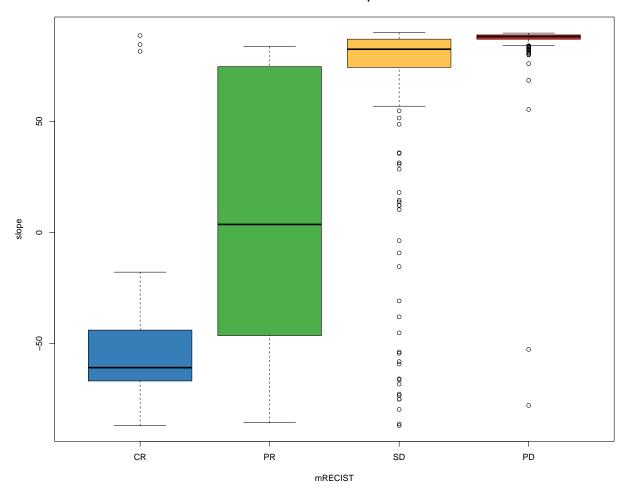

#dev.off()